



QUFORA® URINE BAGS

URINE METERS AND ACCESSORIES

- ➔ Leg bags
- ➔ Bed bags
- ➔ Urine Meters
- ➔ Catheter Valves
- ➔ Accessories





CONTENT


- 3 Qufora A/S
- 4 **Qufora®** Leg bags
- 6 **Qufora®** Bed bags
- 8 **Qufora®** Urine Meters
- 10 **Qufora®** Catheter Valves
- 11 **Qufora®** Accessories

QUFORA A/S

Qufora develops medical devices in the field of gastroenterology and urology. We are a Danish product development company (established in 2006) focusing on innovative medical devices, marketed under the Qufora® trademark.

We are in close dialogue with users and professionals to integrate their knowledge and experience in the product development of new innovative products.

The collaboration has enabled us to develop an exciting portfolio of urine bags and accessories, with focus on product quality and ease-of-use for both professionals and users.

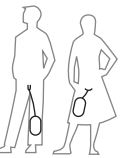

QUFORA® LEG BAGS

Qufora® Leg bags are comfortable to wear and have a nice fit to the leg. The bag is almost invisible underneath the clothes when filled and the soft backing ensures high comfort, even in warm conditions.

Qufora® Leg bags are subject to strict quality requirements to meet the expectations of a user-friendly and durable bag.

PRODUCT FEATURES

- ▶ Five different sizes – 500 ml and 750 ml with a multi-chamber system as well as 350 ml, 600 ml, 750 ml and 800 ml without multi-chamber system.
- ▶ Sterile and non-sterile versions. Most of the bags have an adjustable tube.
- ▶ Sterile Leg bags are available with needle-free outlet for urine sampling in four sizes.
- ▶ The universal connector makes it easy to connect the Leg bag to a catheter or a male external catheter.
- ▶ The discreet, soft and flexible tube ensures a high degree of freedom of movement.
- ▶ The non-return valve prevents the urine from running back in the tubing, thus preventing infection risk.
- ▶ The multi-chamber system ensures that the Leg bag is almost invisible under clothing, even when full.
- ▶ The back side has a skin-friendly ‘non-woven’ soft backing for increased comfort.
- ▶ Available with either the discrete Q-Flow valve, T-valve or Lever valve - all can be operated single-handed.
- ▶ The drainage valves are easy to operate and can be emptied in a hygienic and safe manner.
- ▶ The Leg bag can be connected directly to a Bed bag. This reduces the risk of infection in the event of rupture between catheter and Leg bag.
- ▶ The Leg bags are individually packaged.


Flexible tube with non-return valve



Holes for attaching velcro strap



Scale and comfortable soft backing on multi chamber bags

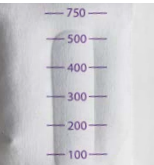

Q-flow or T-tap valve, operated with one hand



Connectors - available in two versions - with or without needle-free sampling



NON STERILE

MULTI CHAMBER - DRAINABLE - ADJUSTABLE TUBE / LOOSE CONNECTOR / CAP - SINGLE PACKED

REF	Volume	Tube	Valve	Soft backing	Qty/box
11121501	500 ml	50 cm	T-valve	Yes	10
11161502	500 ml	50 cm	Q-Flow	Yes	10
11321501	750 ml	50 cm	T-valve	Yes	10
11361502	750 ml	50 cm	Q-Flow	Yes	10

*with 1 pair of straps; ** no adjustable tube.

STANDARD - DRAINABLE - SINGLE PACKED

REF	Volume	Tube	Valve	Soft backing	Qty/box
11061901 *	350 ml	30 cm	Q-Flow	Yes	10
11221502 *	600 ml	50 cm	Q-Flow	Yes	25
11221902	600 ml	30 cm	Q-Flow	Yes	25
11223802 **	600 ml	10 cm	Q-Flow	No	25
11251502 *	600 ml	50 cm	T-valve	Yes	25
11251902	600 ml	30 cm	T-valve	Yes	25/Bag
11259802 **	600 ml	10 cm	T-valve	No	25
11451502	750 ml	50 cm	T-valve	No	25
11451803	750 ml	10 cm	T-valve	No	25
11421502 *	800 ml	50 cm	Q-Flow	Yes	25
11451502 *	800 ml	50 cm	T-valve	Yes	25

*adjustable tube / loose connector / cap; **slanting tube.

COMBI - DRAINABLE - LOOSE CONNECTOR / CAP - SINGLE PACKED

REF	Volume	Tube	Valve	Soft backing	Qty/box
11751101 *	1.500 ml	90 cm	T-valve	Yes	10

*adjustable tube / loose connector / cap

STERILE

DUO – DRAINABLE - NEEDLE-FREE SAMPLE PORT - SINGLE PACKED


REF	Volume	Tube	Valve	Soft backing	Qty/box
12561001 *	1.000 ml	90 cm	Q-Flow	Yes	20

* extra needle-free sample port; “doghouse”; fixation string; hanger; odor filter; wide inlet tube

MULTI CHAMBER – DRAINABLE - NEEDLE-FREE SAMPLE PORT - SINGLE PACKED

REF	Volume	Tube	Valve	Soft backing	Qty/box
12061901 *	350 ml	30 cm	Q-Flow	Yes	10
12161502 *	500 ml	50 cm	Q-Flow	Yes	10
12161802	500 ml	10 cm	Q-Flow	Yes	10
12221502 *	600 ml	50 cm	Lever	Yes	10
12251902 *	600 ml	30 cm	T-valve	Yes	10
12361502 *	750 ml	50 cm	Q-Flow	Yes	10
12361802	750 ml	10 cm	Q-Flow	Yes	10

* adjustable tube / loose connector/cap; 12221502: no needle-free sample port; 12251902 and 12161802: wide inlet tube

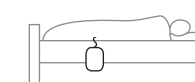

QUFORA® BED BAGS

Qufora® Bed bags are designed for hospitals, nursing homes and home care with focus on functionality and security.

Qufora® Bed bags are available in several versions and sizes to suit individual needs.

PRODUCT FEATURES

- ▶ A wide product assortment. Available in sterile or non- sterile, long or short tubes, drainable and non-drainable and with or without needle-free sampling.
- ▶ Has a universal connector which easily connects to the catheter or Leg bag.
- ▶ The flexible connection tube provides for safe drainage and allows the bag to be placed below bladder level.
- ▶ The non-return valve prevents the urine from running back in the tubing, thus preventing infection risk.
- ▶ Capacity: 750ml, 1.500 ml or 2.000 ml.
- ▶ Two different drainage valves - free-flow T-valves or large Lever valve. The drainage valves are easy to operate.
- ▶ The Bed bag can be emptied in a hygienic and safe way.
- ▶ Clear scale on the outside makes it easy to read the urine volume.
- ▶ The free-flow or large outlet ensures that clots freely can pass.
- ▶ Tube lengths from 10 cm - 130 cm.

“Eyes” for attaching to a bed bag stand



Connectors - available in two versions - with or without needle-free sampling



Both drainage valves are easy to operate with one hand



Non-drainable bag



“Doghouse”



NON STERILE

NON-DRAINABLE – TEAR MARK – SINGLE PACKED

REF	Volume	Tube	Valve	Single packed	Qty/box
13303802	750 ml	10 cm	-	Yes	25
13703412	1.500 ml	110 cm	-	Yes	25
13703802	1.500 ml	10 cm	-	Yes	25
13703902	1.500 ml	35 cm	-	Yes	25
13800101	2.000 ml	90 cm	-	No	10
13803112	2.000 ml	90 cm	-	Yes	25

DRAINABLE – SINGLE PACKED

REF	Volume	Tube	Valve	Single packed	Qty/box
13753101	1.500 ml	90 cm	T-valve	Yes	25
13850301	2.000 ml	130 cm	T-valve	Yes	10
13853101	2.000 ml	90 cm	T-valve	Yes	25
13853401	2.000 ml	110 cm	T-valve	Yes	25
13853402 **	2.000 ml	110 cm	T-valve	Yes	100
13853403 **	2.000 ml	110 cm	T-valve	Yes	100

** wide inlet tube

STERILE

DRAINABLE – NEEDLE-FREE SAMPLE PORT – ODOR FILTER - SINGLE PACKED

REF	Volume	Tube	Valve	Single packed	Qty/box
14823202 – U4	2.000 ml	120 cm	Lever	Yes	25
14851201 – S4	2.000 ml	120 cm	Lever	Yes	25
14853202 – Q4	2.000 ml	120 cm	T-valve	Yes	25
16823201 *	2.000 ml	130 cm	Lever	Yes	20

- U4/Q4: fold-away holder; bed sheet clip; tube clamp; wide inlet tube

- S4: “dog-house”; bed sheet clip; tube clamp; wide inlet tube

* closed system: drip chamber; “doghouse”; hook; odor filter; fixation string; bed sheet clip

QUFORA® URINE METER CLOSED SYSTEM

The sterile Qufora® Urine Meter Closed System Compact is a specially designed urine collection bag system for professional use. The Urine Meter is mainly used in intensive care units and post-operative wards where urine output is being closely monitored.

PRODUCT FEATURES

- ▶ The Qufora® Urine Meter Closed System Compact comes with a fixed and easy to clean sampling port which allows for needle-free sampling and eliminates the risk of needle injuries.
- ▶ 500 ml uri-flow measuring chamber with Pasteur drip chamber and anti-reflux valve to prevent retrograde migration of bacteria and venting air inlet filter for pressure balance.
- ▶ Small volumes of urine can be monitored - useful for pediatric use.
- ▶ Built-in air valve system in the flow chamber provides permanent ventilation.
- ▶ Integrated bacteria-proof and liquid-tight ventilation filter ensures an optimal pressure balance.
- ▶ The kink-resistant inlet tube ensures a reliable drainage of urine from the bladder into the measuring chamber.
- ▶ The translucent films are extra thick for added durability.
- ▶ The Doghouse fixates the outlet tube (fold-away position) to prevent contamination and to keep distance to the floor. A simple one hand operation.
- ▶ The fixation string (75 cm long) and integrated hanger enable fastening to the bed frame.
- ▶ The sheet clip is provided for easy fixation of the tube to the bed linen to reduce the risk of the tube kinking.
- ▶ The sliding clamp is useful when taking a urine sample or changing the bag.
- ▶ The protection cap prevents urine from dropping from the inlet tube during exchange of bags.
- ▶ The Lever valve has a free-flow outlet that ensures that clots can freely pass.
- ▶ The flexible connection tube provide for safe drainage and allows the bag to be placed below the bladder level.



Connector with needle-free sample-port



Doghouse



STERILE

URINE METER – REPLACEMENT URINE BAG

REF	Volume	Tube length	Quantity poly bag	Quantity shipper box
16923201	2.600 ml bag + 500 ml chamber	130 cm	10	10

QUFORA® URINE METER MEASURING SYSTEM

The sterile Qufora® Urine Meter Measuring System is a specially designed urine collection measuring system for professional use. The Urine Meter Measuring System is mainly used in intensive care units and post-operative wards where urine output is being closely monitored. The system has built-in safety features and easy to use functionalities.

PRODUCT FEATURES

- ▶ The Qufora® 500 ml Urine Meter Measuring System with three (3) chambers has a precise and easy-to-read scale in graduations of 1 ml to 5 ml to 10 ml.
- ▶ The built-in safety features include the vented drip chamber with the integrated non-reflux valve, overflow chamber, anti-kink tube above the chamber and the bayonet locking mechanism for attaching the bag.
- ▶ Features include a built-in liquid-tight air ventilation filter in the chamber, a needle-free sampling port with a protection cap, strings for hanging the chamber and a tube clamp and bed sheet clip.
- ▶ The 2.000 ml bag includes a non-return valve, free-flow outlet T-valve, hydrophobic air vent filter and the valve fold-away.



NON STERILE

URINE METER – REPLACEMENT URINE BAG

REF	Volume	Measurement Scale	Free-Flow outlet valve	Non-return valve	Hydrophobic Air filter	Quantity Retail / Shipper unit
19853201	Bag: 2.000 ml	100 – 2.000 ml: 100 ml	Yes	Yes	Yes	50 Non-Sterile Replacement Bags

STERILE

URINE METER – URINE MEASURING SYSTEM – SINGLE PACKED

REF	Volume	Tube length	Measurement Scale	Needle-Free Sample Port	Chamber Lock	Quantity Retail / Shipper unit
18853201	Chamber: 500 ml	150 cm	5 – 50 ml: 1 ml 60 – 120 ml: 5 ml 130 – 500 ml: 10 ml	Yes	Yes – Bayonet system for quick exchange of the bag	10 Sterile Systems
-	Bag: 2.000 ml	-	100 – 2.000 ml: 100 ml	-	-	

QUFORA® CATHETER VALVES

The catheter valves can be attached to permanent catheters as an alternative to a urine bag, which provides the user more freedom and discretion.

All valves are easily operated with one hand.

The valves can be connected to Leg and Bed bags.

Single packed. Sterile.

PRODUCT FEATURES

- ▶ The Q-Flow valve has a soft silicone tube and is easy to put together with a urine bag or catheter.
- ▶ The Q-Flow valve is easy to handle with one hand - also for people with reduced hand mobility. The valve secures a free flow of urine at between 90° and 180°. Furthermore, 180° ensures that the valve does not close during the night.



- ▶ T-valve and Lever valve, both equipped with a soft silicone tube and easy to handle with one hand.

Q-Flow valve



T-valve



Lever valve



QUFORA® VALVES

REF	Description	Qty/box
21104101	T-valve	30
21104201	Q-Flow valve	30
21104301	Lever valve	30

QUFORA® ACCESSORIES

Qufora® Accessories are designed for easy and discreet use to give the user greater freedom of movement.

PRODUCT FEATURES

- ▶ Hanger. It is durable and fits all beds and can easily be operated with one hand.
- ▶ Velcro strap for attaching the Leg bag to the leg, so it can be worn discreetly under clothing.
- ▶ The floor stand can be used for both Leg bag and Bed bags. On one side, the floor stand is equipped with one large hook for a water bag and the other side with four small hooks for a bed bag.

Hanger



Velcro strap



Floor stand



QUFORA® ACCESSORIES

REF	Description	Qty/box
21100101	Bed hanger	5
21100301	Leg straps – 1 pair	10
99910009	Floor stand	1



At Qufora A/S we know what challenges our users face in everyday life. Thanks to our extensive experience and in-depth product knowledge, we know of the hidden problems and complications that many people think they are the only ones having.

Based on this knowledge, our vision is to create products that help provide healthy lives for people with chronic conditions by breaking the limits of conventional thinking.

For us, it is about adding quality of life and facilitating everyday life for people - both for the caregiver in their daily work in hospitals and other institutions and for the user in his or her everyday life.

CONTACT

Qufora A/S
Gydevang 28-30
DK-3450 Allerød
Denmark

Tel: +45 88 30 05 92
info@qufora.com
www.qufora.com